ALE-HOP

User manual | Manual de usuario | Manual
d'usuari | Manual do usuário | Mode
d'emploi | Benutzerhandbuch | Manuale
d'uso | Priručnik s uputama | للتحويل بين الأوضاع

MASSAGE GUN HOT AND COLD

Product specifications:

Voltage: 7.4V

Battery capacity: 1800 mAh

Rated power: 4.7W (Normal) 7.2W (Heat)

Charging Time: +/- 3-4 Hours

Working Time: +/- 1-1.5 Hours

Working Current (Max.): 500mA


(Normal) 800mA (Heat)

Heating Temperature: 60°C (±5°C)
in 3 minutes

Rated time: 15 minutes

Speed at Mode 1: 2800-2900rpm	Speed at Mode 2: 2900-3000rpm
Speed at Mode 3: 3000-3200rpm	Speed at Mode 4: 3200-3400rpm
Speed at Mode 5: 3400-3600rpm	Speed at Mode 6: 3600-3800rpm

Product overview:

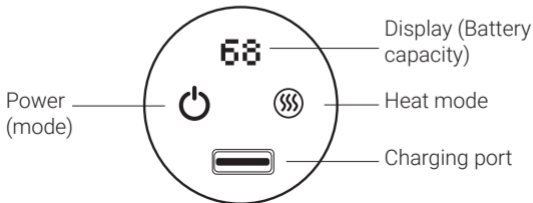Selected mode
display

Heat mode is not
activated

Heat mode is
activated

Hot head	Cold head	Round head
 <p data-bbox="91 386 324 458">To soothes stiff, aching, muscles</p>	 <p data-bbox="399 386 640 495">For joint swelling and muscle inflammation</p>	 <p data-bbox="710 386 951 458">For large muscle groups</p>
Flat head	Fork head	Bullet head
 <p data-bbox="94 873 324 945">For any parts of body</p>	 <p data-bbox="401 873 638 904">For neck & spine</p>	 <p data-bbox="762 873 899 904">For joints</p>

Main features:

- Rechargeable cordless massager.
- Massage mode with 6 power levels.
- Hot and cold therapy.
- Helps relieve muscle stiffness and soreness.
- Improves range of motion.
- Promotes blood circulation.
- Powerful high torque motor.
- Lightweight and easy to use for myofascial self-release.

Operation:

1. Choose a suitable massage head and plug it into the massage gun.
2. Press and hold the power button (Mode) for 1-2 seconds to turn it on.
3. Press briefly to choose another massage mode.
4. Switches to 6 different massage modes.
5. Automatic switch-off after approx. 15 minutes.
6. Adjust the massager to the desired speed on your body and apply the appropriate pressure.
7. When massage movement starts to slow down, stop using the massager, and recharge it by using the supplied

charger cable.

8. Press and hold the Power (Mode) button for 1-2 seconds to turn off the unit.

Note: Press the Power (Mode) button to display the remaining battery capacity when the product is switched off.

How to use Hot Head:

The Hot Head attachment helps relieve muscle tension and pain.

Plug the Hot Head attachment into the massage gun, press and hold the power (mode) button for 1-2 seconds to turn the unit on, then press and hold the heat mode button for 1-2 seconds to activate the heat mode. The red light will illuminate on the heat mode button.

Wait 3-4 minutes for the head to heat up. Position the massage gun at the desired speed on the desired spot.

Remember to turn off the heat mode by long pressing the heat mode button for 1-2 seconds after the massage.

How to use Cold Head:

The cold head helps to treat muscle inflammation. Before using the cold head, place it in the freezer for about 2 hours for best results.

Connect the cold head to the massager and turn it on (press and hold the power (mode) button for 1-2 seconds). Adjust the massage gun to the desired speed and at the desired location.

How to charge:

- Before using it for the first time, charge it fully (about 3-4 hours).
- Connect the USB charging cable (included) to a suitable USB socket with DC power.
- The battery capacity indicator will light up and go up when the battery is charging.

How to clean:

- Before cleaning, make sure that the power switch is turned off and the power plug is disconnected from the product.
- Use a soft cloth to wipe the surface of the product and the control panel.

- Use a slightly damp cloth for cleaning and then wipe it dry.

Safety precautions:

- For adults only
- Use only on the soft tissues of the body to avoid the pain or discomfort. Do not use it on head, bones, or other hard areas
- Keep fingers, hair or other body parts from the shaft and the back of the massager head to avoid pinching.
- Do not place anything into the ventilation holes on the shaft.
- Do not immerse this product in water or allow water to enter the product.
- Do not drop or misuse this product.
- Keep children and pets away from the movable parts of this product.
- Do not use attachments not recommended.
- Do not use this product outdoors.
- Do not use when the ambient temperature is above 40°C and below 5°C.
- The recommended usage time of this product is 15 minutes/session.
- Check if any part of this product is damaged before use. If so, do not use it.
- Do not fall asleep while using this

product.

- Use only on a dry surface. Lightly press and move over the desired areas, massaging for approximately 60 seconds on each area.
- This product is not intended to be used for medical treatment purposes.
- At the end of its useful life, dispose it in the proper container.

People with the following symptoms should consult a doctor before use:

- People suffering from skin conditions such as inflamed skin, sunburn, or other dermatology diseases.
- People suffering from any blood circulation disorders such as thrombosis, severe aneurysms, or other skins infections.
- People suffering from osteoporosis.
- People suffering from heart disease or wearing electronic medical equipment such as pacemakers.
- Pregnant women.

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DÉNIA, S.A.U. | A-03740966 | Avda. Alacant, 174, Bellreguard.