# ALE-HOP

**User manual** | Manual de usuario | Manual  
d'usuari | Manual do usuário | Mode  
d'emploi | Benutzerhandbuch | Manuale  
d'uso | Priručnik s uputama | للتحويل بين الأوضاع

**EN**

# **MICROCURRENT FACIAL MASSAGER**

## **PRODUCT SPECIFICATIONS**

Rated power: 0.15W

Output voltage: 5V

Working voltage: 2.4V

Battery capacity: 800mAh


Use time: 2-3 hours

Charging time: 2 hours

## **CHARGING THE INSTRUMENT**

Please charge the device for 12 hours before its first use. To charge the device, connect one end of the provided charging cable to a USB power source and the other to the device. When the device is charging, the LED charging indicator lights will increasingly illuminate. When it is fully charged, the indicator lights will stay illuminated. The device will not turn on while it is charging.

## **COMPONENT DESCRIPTION**

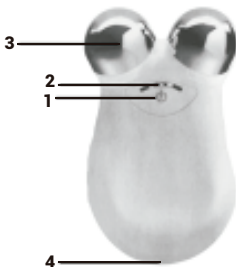1. Power ON/OFF and Intensity control UP/DOWN button
2. Intensity level indicator
3. Micro current treatment electrodes
4. USB interface

## USE STEPS

1. Please check the "Contraindications of use", "Precautions", and "Warnings" sections before using the device.
2. Please clean your face with an oil-free cleanser. Men need to shave before care because facial hair may interfere with the conductivity of the spheres.
3. Please apply gelatine or emulsion to clean and dry skin. Always use gel or emulsion when using the facial massager device.
4. Please adjust the intensity level by gradually pressing the ON/OFF button to increase or decrease the microcurrent output. When you switch among the five intensity levels, the device will beep, and the number of illuminated blue lights will change.


### Intensity levels:

**Minimum intensity:** left light on / left two lights on

**Medium intensity:** middle light on

**Maximum intensity:** right light on / right two lights on

5. Place the device directly in contact with the treatment area. Always keep both spheres in contact with the skin to keep the microcurrent flowing. Using light to medium pressure, begin the treatment. Facial skin sensitivity varies from person to person. Do not worry if you do not feel anything on your skin. If you experience tingling, add more gel or lotion to the area you are treating. Note: the facial massager will not contract muscles. When used correctly, the current is gentle and very soothing.

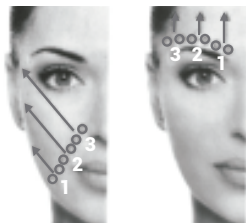6. After the treatment is completed, wipe the gel or lotion from the device with a damp cloth and massage the remaining gel or lotion into the skin.

### **TREATMENT FREQUENCY**

- Use 5 minutes a day.
- For the first 60 days, use 5 days per week.
- After 60 days, use 2-3 times per week for maintenance.
- Allow 24 hours between treatments.
- Do not use more than twice a day.

### **MAINTENANCE AND SERVICE**

1. Please turn off the instrument before cleaning.
2. Use a slightly damp, soft cloth with water or rubbing alcohol to clean the device after every treatment. Protect the instrument from moisture. When cleaning, do not immerse the device in water to prevent water from entering the charging port. Do not put the instrument into a dishwasher or microwave oven.
3. Do not use sharp or corrosive, volatile diluent cleaners.

### **CONTRAINDICATIONS FOR USE**

- Do not use this device if you are pregnant. The safety of electrical stimulation during pregnancy has not been established.
- Do not use on infants and minors under the age of 18.
- Do not use it if you have an electronic implanted device such as a defibrillator, neurostimulator, pacemaker, or ECG Monitor. If you use a body-worn medical device such as an insulin pump or other electronic device, consult a healthcare professional before using this device.
- Do not use if you have any active cancer or any suspicious or cancerous lesions, especially skin precancerous nevus (carnivorous nevus)

- Do not use this device if you have a history of epilepsy or seizures, heart disease or immune system damage.
- Do not use it if you have a skin disease such as sores, psoriasis, eczema, rash, recurrent herpes simplex, etc.
- Apply stimulation only to normal, intact, clean, healthy skin. Do not use over swollen, infected, or inflamed areas or skin eruptions, e.g., phlebitis, thrombophlebitis, broken capillaries, varicose veins, scars, wounds, sensitive skin etc.
- Please avoid using the device after undergoing a laser skin rejuvenation process, chemical peeling, or wound healing. Use the device after the skin is completely healed. If there are more medical problems, please consult your healthcare provider.
- Please do not use this device if you have been implanted recently temporary fillers (hyaluronic acid, water light needle, protein thread) and permanent fillers (prosthesis nose augmentation and chin augmentation). Only use it after more than three months after the intervention and after consultation with your doctor.
- Do not use if you have a metal allergy.
- In the process of using, please stop using immediately in case of discomfort.

## **PRECAUTIONS**

- The long-term effects of transcutaneous electrical stimulators are unknown.
- If you have suspected or diagnosed heart disease, you should follow precautions recommended by your healthcare provider.
- Some users may experience skin irritation or hypersensitivity due to the electrical stimulation. Irritation can usually be reduced by using a gelatine or emulsion with the device.
- Skin irritation and burns beneath the electrodes have

been reported with the use of other powered muscle stimulators. If you experience skin irritation or burning, discontinue use immediately.

- If using the facial massager device causes persistent pain or discomfort, discontinue use immediately.
- A slight tingling sensation may occur while using the device. Lowering the intensity may reduce or eliminate this sensation.
- There may be a perception of flashing lights during the use of the facial massager, resulting from the stimulation of the eye area. If you notice this condition constantly without the use of your device, consult your healthcare provider.
- Some users may experience a metallic taste in their mouth during the use of your device, resulting from stimulation and varying sensitivity around the mouth. If you notice this condition constantly without the use of your device, consult your healthcare provider.
- The battery of the device is not removable by the user. If the device is to be disposed of, do not put into a fire. Dispose of the device per local regulations.
- Inspect the device and do not use the device if you see any signs of internal battery leaks.
- Store and use the device only in a room temperature environment. Do not store it in the sun or on a hot surface. High temperatures may be damaging. The device should be operated, transported, and stored in the manner specified in this manual. If the instrument is abnormally hot, please stop using and contact the factory for after-sales service.
- Only use the power adapter provided with this device.

## **WARNINGS**

Please read this manual carefully before use and save it for future reference.

1. Please follow the operation instructions. Do not use

the machine for other purposes.

2. The device is designed for cosmetic and personal use. If this is not the case, please consult a healthcare professional before use.
3. Before using the device, consult a healthcare professional if you have any underlying skin conditions.
4. The product may not be safe for use on the neck in individuals with a known history of carotid artery disease, stroke, transient ischemic attack (TIA), carotid stenosis, unstable blood pressure, easy fainting, thyroid disease, or prior surgical procedures in the neck area.
5. Do not use it in the eye area (circular muscle within the orbital rim), on the temples, the breast or the groin area.
6. Do not apply stimulation across your chest. The introduction of electrical current into the chest may cause rhythmic disturbances to your heart, which could be lethal.
7. Do not use while driving, operating machinery or during any activity in which involuntary muscle contractions may put the user or the public at unnecessary risk of injury.
8. Never operate the device if it's not functioning as always or not operating properly due to damage caused by dropping or accidental submersion in water. If the instrument accidentally falls into the water, please shut it down immediately. Do not try to start it and contact the factory after-sales service.
9. Never drop or insert any object into any openings on the device.
10. The performance of procedures other than those specified in the instruction manual may result in injury.
11. Avoid applying excessive pressure when using the device on the anterior neck.
12. If you experience discomfort, adverse reaction or difficulty breathing or swallowing from using the device on the neck, discontinue use.

13. If you experience symptoms of hypotension (dizziness or light-headedness) while using the device, discontinue use.
14. Keep all electrical appliances away from water (including baths, showers, and sinks). Do not apply stimulation while bathing or showering.
15. Care should be taken to minimize treatment directly over the midline of the neck, specifically avoiding the thyroid.
16. If the device begins to overheat, turn the device OFF immediately.
17. Do not apply stimulation over painful areas. Consult a healthcare professional before using this device.
18. Supervise children to ensure they do not play with the instrument.
19. Put the packaging materials (plastic bags) in places children cannot reach.
20. If the instrument is damaged, to avoid danger, it must be repaired by the manufacturer and professional maintenance personnel.
21. Do not drop or knock the instrument. If the instrument is damaged due to dropping or collision, do not use it.
22. Because the instrument is equipped with a polymer battery, do not place the instrument in a high-temperature environment or the sun.
23. Do not use cables other than those specified as replacement parts for internal components. It may result in increased emissions or decreased immunity of the system. Do not use near the thyroid gland.
24. Do not operate the device in the presence of electronic monitoring equipment (e.g., heart monitors, heart tracing, or in proximity (1 meter) to shortwave, microwave, or electromagnetic therapy equipment, as it may cause the medical equipment to malfunction. Such equipment may produce instability in the device and can cause burns at the stimulation site.
25. Do not use this device when there is a tendency to bleed following acute injury or breaks in bones.

26. Do not apply stimulation over areas of skin that lack normal sensation.
27. The device should not be used at extreme temperatures.
28. Keep this device out of the reach of children and pets.
29. Discontinue using the device immediately if you suspect the device to be malfunctioning, or not operating properly.
30. Do not disassemble the device as this may cause damage, malfunction, electrical shock, fire, or injury.

### **BOX CONTENT**

- 1x Facial massage device
- 1x Charging cable
- 1x Clothing bag
- 1x Cleaning cloth
- 1x Instruction manual

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**DÉNIA, S.A.U.** | A-03740966 | Avda. Alacant, 174, Bellreguard.